

# Installation Instructions for Radiant Heater RH600

**Dear Customer,**

Before shipping we have completed a Periodic Manufacturer's Service including the change of all parts for a 4 year PMS on your device. The following parts are renewed:

- Battery P/N 83 01 856
- Time Keeper RAM P/N 18 28 142

Please install the Radiant Heater RH600 with the delivered screws and washers according to the following sections of this document. If you have any questions, please feel free to call Dräger Medical, Inc. at 1 800 4 DRAGER and ask for Service.

Sincerely,

Andreas Lenke  
Technical Product Manager  
DrägerService

## 1 Installation of RH600 on the Dräger wall or stand mount




**The RH600 is heavy, please ask a second person for assistance.**

- Install the RH600 by using the 4 screws M4x12 DIN 912 (P/N 12 93 664) with the 4 spring washer B4 DIN 127 (P/N 13 31 345) for this installation, do not use the 4 washers A4.3 DIN 125 (P/N 13 27 542). Secure the 4 screws with Loctite 221 or 222.



**Secure the 4 fixing screws with Loctite 221 or Loctite 222. For your convenience we have sent one bottle of Loctite with the Radiant Heater.**

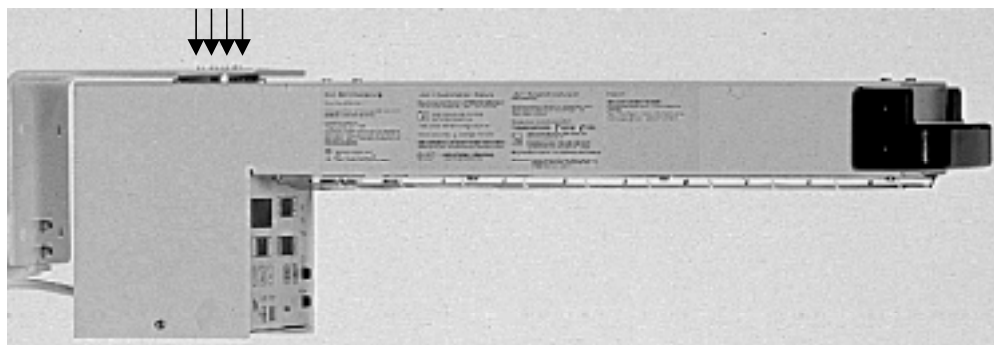


## 2 Installation of RH600 on a Babytherm 8000 OC:



The RH600 is heavy, please ask a second person for assistance.


- Unscrew 3 screws and remove cover plate from the Babytherm 8000 OC.

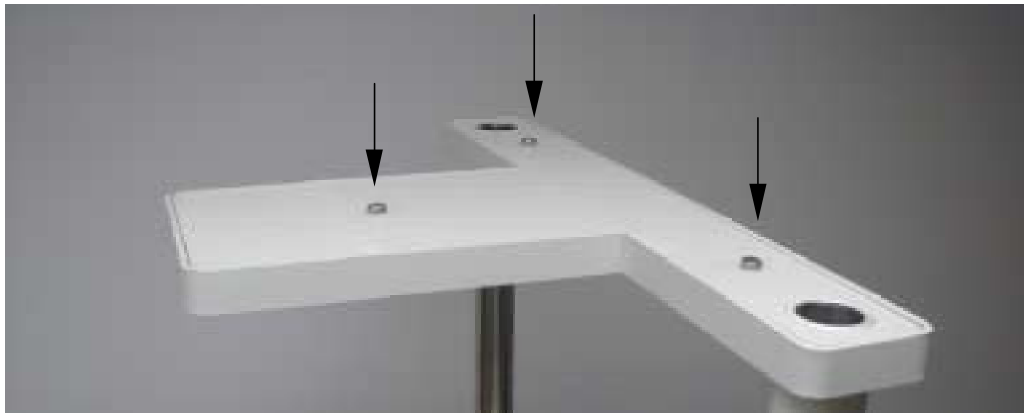

- Install the RH600 by using the 4 screws M4x12 DIN 912 (P/N 12 93 664) with the 4 washers A4.3 DIN 125 (P/N 13 27 542) and the 4 spring washer B4 DIN 127 (P/N 13 31 345). Secure the 4 screws with Loctite 221 or 222.



**Secure the 4 fixing screws with Loctite 221 or Loctite 222. For your convenience we have sent one bottle of Loctite with the Radiant Heater.**

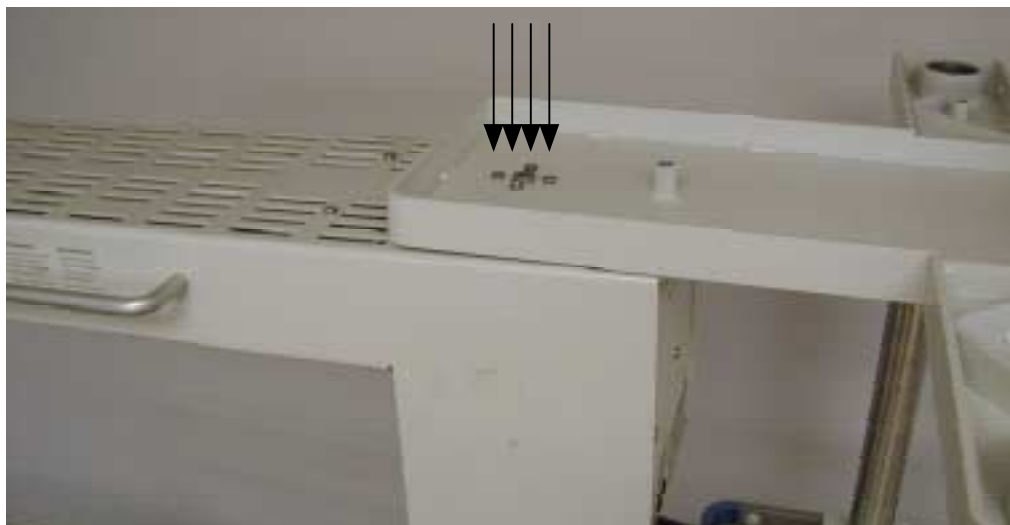

- Assemble cover plate with the 3 screws.



End of Instructions.